

PAB215Ra01**Polyclonal Antibody to Fibrinogen Beta (FGb)****Organism Species: Rattus norvegicus (Rat)****Instruction manual**

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PRODUCT INFORMATION]**Immunogen:** FGb, Rat**Clonality:** Polyclonal**Host:** Rabbit**Immunoglobulin Type:** IgG**Purification:** Affinity Chromatography.**Applications:** WB, ICC, IHC-P, IHC-F, ELISA**Concentration:** 200µg/mL**UOM:** 100µg**[IMMUNOGEN INFORMATION]****Immunogen:** Recombinant FGb (Gly33~Gln479) expressed in *E.coli*.**USCN Accession No.:** RPB215Ra01**Sequence:** The target protein is fused with two N-terminal Tags, His-tag and S-tag and its sequence is listed below.

MHHHHHSSG LVPRGSGMKE TAAAKFERQH MDSPDLGTDD DDKAMADIGS EF-
GHRPVDRR KEEPPSLRPA PPPISGGGYR ARPAKVDAGQ KKVERKPPDA GGCVHGDGDM
GVLCP TGCEL RQTLNHERP IKNSIAELNS NINSVSETSS VTFQYLTLK DMWKKKQAQV
KDNENVINEY SSILEDQKLY IDETVNDNIP LNLRLRSIL EDLRSKIQKL ESDISAQTEY
CHTPCTVNCN IPVVS GKECE EIIRKGGETS EMYLIQPDTS SKPYRVYCDM KTEGGWTVI
QNRQDGSVDF GRKWDPYKKG FGNIATNEDT KKYCGLPGEY WLGNDKISQL TRIGPTTELLI
EMEDWKGDKV KAHYGGFTVQ TEANKYQVSV NKYKGTAGNA LMEGASQLVG
ENRTMTIHNG MFFSTYDRDN DGWVTTDPRK QCSKEDGGGW WYNRCHAANP
NGRYYWGGLY SWDMSKHGTD DGVVWMNWKG SWYSMRRMSM KIRPVFPQQ

[ANTIBODY SPECIFICITY]

The antibody is a rabbit polyclonal antibody raised against FGb. It has been selected for its ability to recognize FGb in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500


Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200


Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

[IMAGES]

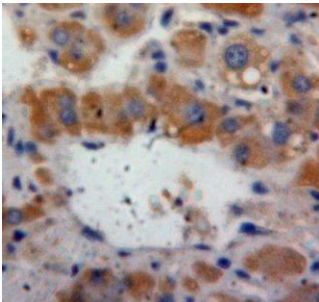

Used in DAB staining on formalin fixed
paraffin-embedded Liver tissue

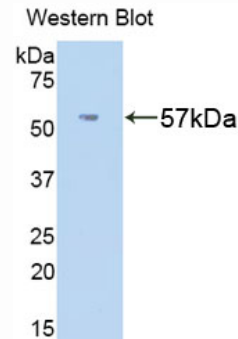

Used in Western Blot, Sample:
Recombinant FGb, Rat

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN_3 , 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant FGb (Gly33~Gln479) dispensed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.
5uL per well when used in enhanced chemiluminescent (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN_3 0.02%.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.