



A90775Ra01 Polyclonal Antibody to Hexokinase 1 (HK1) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

5th Edition (Revised in January, 2013)

[PRODUCT INFORMATION]

Immunogen: HK1
Clonality: Polyclonal


Host: Rabbit

Species Reactivity: Rat

Ig type: Rabbit IgG

Purification: Antibodies are purified by target protein affinity chromatography. **Applications:** WB, ICC, IHC-P, IHC-F

Form: Liquid Size: 100µg

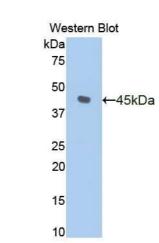

Sample: Recombinant rat HK1

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant rat HK1 (Asp413~Val540) expressed in E.coli.

Molecular Weight: 41.7kDa

USCN accession No.: P90775Ra01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and

GST-tag and its sequence is listed below.

MRNKKFELGL EFPNLPYYID GDVKLTQSMA IIRYIADKHN MLGGCPKERA EISMLEGAVL DIRYGVSRIA YSKDFETLKV DFLSKLPEML KMFEDRLCHK TYLNGDHVTH PDFMLYDALD VVLYMDPMCL DAFPKLVCFK KRIEAIPQID KYLKSSKYIA WPLQGWQATF GGGDHPPKSD GSTSGSGHHH HHHSAGLVPR GSTAIGMKET AAAKFERQHM DSPDLGTLEV LFQ GPLGSEF- DGSLYKMH PQYSRRFHKT LRRLVPDSDV RFLLSESGTG KGAAMVTAVA

11271 Richmond Avenue, H104, Houston, TX 77082, USA | Toll free: 001-888-960-7402 | Fax: 001-832-538-0088 | Email: mail@uscnk.us | Http://www.uscnk.us Export Processing Zone, Wuhan 430056, P.R. China | Toll free: 0086-800-880-0687 | Fax: 0086-27-8425-9551 | Email: mail@uscnk.com | Http://www.uscnk.com





YRLAEQHRQI EETLAHFRLS KQTLMEVKKR LRTEMEMGLR KETNSKATVK MLPSFVRSIP DGTEHGDFLA LDLGGTNFRV

[ANTIBODY SPEFICITY]

Anti HK1 is a rabbit polyclonal antibody raised against rat HK1. This antibody has been selected for its ability to recognize rat HK1 in immunohistochemical staining and western blotting, non cross-reactive with other members of the family.

[APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200

Optimal working dilutions must be determined by end user.

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH 7.4, containing 0.02%Na₃N, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant rat HK1 (Asp413~Val540) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading buffer: 100mM Tris(PH 8.8), 2%SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, Na₃N 0.02%.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.