

**EPD025Ge51 1mg**  
**Eukaryotic Green Fluorescent Protein (GFP)**  
**Organism Species: *Pan-species (General)***  
***Instruction manual***

FOR RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

12th Edition (Revised in Aug, 2016)

**[ PROPERTIES ]**

**Source:** Eukaryotic expression

**Host:** Yeast

**Residues:** Met1~Lys238

**Tags:** N-terminal His Tag

**Subcellular Location:** Secreted

**Purity:** > 95%

**Traits:** Liquid

**Buffer formulation:** PBS, pH7.4, containing 20%Glycerine.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.7

**Predicted Molecular Mass:** 27.1kDa

**Accurate Molecular Mass:** 28kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]****[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ SEQUENCE ]**

MSKGEELFTG VVPILVELDG DVNGHKFSVS GEGEDATYG KLTLKFICTT  
GKLPVPWPTL VTTFSYGVQC FSRYPDHMKQ HFFKSAMPE GYVQERTIFF  
KDDGNYKTRA EVKFEGDTLV NRIELKGIDF KEDGNILGHK LEYNYNSHNV  
YIMADKQKNG IKVNFKIRHN IEDGSVQLAD HYQNTPIGD GPVLLPDNHY  
LSTQSALSKD PNEKRDHMLV LEFVTAAGIT HGMDELYK

**[ IDENTIFICATION ]**

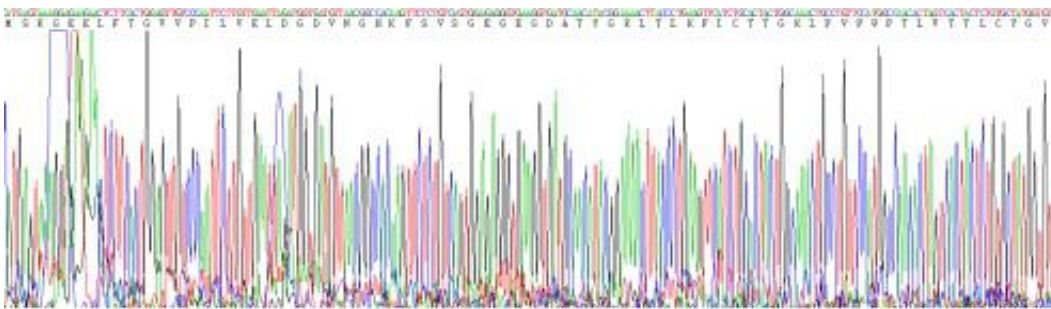

**Figure . Gene Sequencing (extract)**

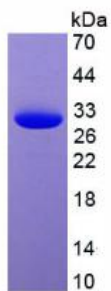


Figure. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.