

#### MAA133Ca21

**Monoclonal Antibody to Tumor Necrosis Factor Alpha (TNFa)** Organism Species: Canis familiaris; Canine (Dog) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

#### [ PRODUCT INFORMATION ]

Immunogen: TNFa, Canine

**Clonality:** Monoclonal

Host: Mouse


Immunoglobulin Type: IgG

**Purification:** Affinity Chromatography.

Applications: WB, ICC, IHC-P, IHC-F, ELISA

Concentration: 500µg/mL

**UOM**: 200µg

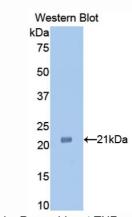

Sample: Recombinant TNFa, Canine

### [ IMMUNOGEN INFORMATION ]

Immunogen: Recombinant TNFa (Leu63~Leu233) expressed in E.coli.

USCN Accession No.: RPA133Ca01

Sequence: The target protein is fused with N-terminal His-Tag and its sequence

is listed below.

MGHHHHHHSGSEF-LPNGLQLI SPLAQTVKSS SRTPSDKPVA HVVANPEAEG QLQWLSRRAN ALLANGVELT DNQLIVPSDG LYLIYSQVLF KGQGCPSTHV LLTHTISRFA VSYQTKVNLL SAIKSPCQRE TPEGTEAKPW YEPIYLGGVF QLEKGDRLSA EINLPNYLDF **AESGQVYFGI IAL** 



#### [ANTIBODY SPECIFITY]

The antibody is a mouse monoclonal antibody raised against TNFa. It has been selected for its ability to recognize TNFa in immunohistochemical staining and western blotting.

# [APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

## [CONTENTS]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

#### [STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.