

#### CPB301Ge21 100µg OVA Conjugated Asymmetrical Dimethylarginine (ADMA) Instruction manual


FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

11th Edition (Revised in May, 2016)

### [PROPERTIES]

Source: Protein Conjugation

**Original Structure:** 

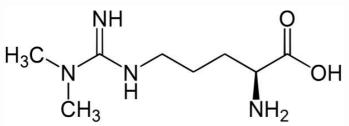

Original Chemical Formula: C<sub>8</sub>H<sub>18</sub>N<sub>4</sub>O<sub>2</sub>

Original Mol. Mass: 202.3Da

**Purity: >90%** 

Traits: Freeze-dried powder.

Buffer Formulation: PBS, pH7.4.

Applications: Immunogen; Coating Antigen; ELISA; SDS-PAGE.

## [ <u>USAGE</u> ]

Reconstitute in PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

# Coud-Clone Corp.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

#### [IDENTIFICATION]

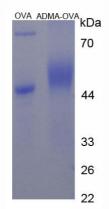


Figure 1. SDS-PAGE

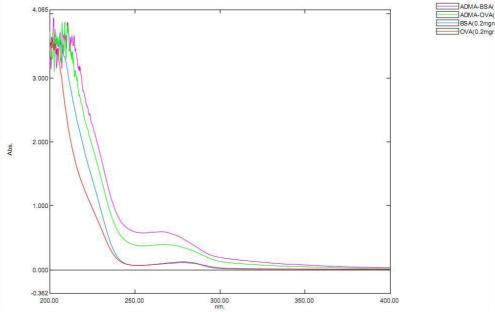


Figure 2. UV-spectrum

1304 Langham Creek Dr. Soire 226, Houston, TX 77084, USA | 001-888-960-7402 | www.cloud-chore.ox | mail@cloud-clore.ox Export Processing Zone, Wuhan, Hubei 430056, PRC | 0086-000-880-0687 | www.cloud-clore.com | mail@cloud-clore.com