PAA704Hu01 Polyclonal Antibody to Apolipoprotein E (APOE) Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PRODUCT INFORMATION]

Immunogen: APOE, Human Clonality: Polyclonal Host: Rabbit Immunoglobulin Type: IgG Purification: Affinity Chromatography. Applications: WB, ICC, IHC-P, IHC-F, ELISA Concentration: 200µg/mL UOM: 100µg

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant APOE (Lys19~Gln181) expressed in E.coli.

Accession No.: RPA704Hu01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and T7-tag and its sequence is listed below.

MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF- KV EQAVETEPEP ELRQQTEWQS GQRWELALGR FWDYLRWVQT LSEQVQEELL SSQVTQELRA LMDETMKELK AYKSELEEQL TPVAEETRAR LSKELQAAQA RLGADMEDVC GRLVQYRGEV QAMLGQSTEE LRVRLASHLR KLRKRLLRDA DDLQKRLAVY Q

[ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against APOE. It has been selected for its ability to recognize APOE in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400


Immunocytochemistry in formalin fixed cells: 1:100-500 Immunohistochemistry in formalin fixed frozen section: 1:100-500 Immunohistochemistry in paraffin section: 1:50-200

Designed by Cloud-Clone Corp., Assembled by Uscn Life Science Inc. ISO9001:2008; ISO13485:2003 11271 Richmond Avenue Suite H104. Houston, TX 77082, USA | Toll free: 001-888-960-7402 | Fax: 001-832-538-0068 | Http://www.cloud-clone.us | E-mail: mail@cloud-clone.us Export Processing Zone Building F. Wuhan, Hubei 430056, PRC | Toll free: 0086-800-880-0667 | Fax: 0086-27-8425-9551 | Http://www.uscnk.com | E-mail: mail@cloud-clone.us

Cloud-Clone Corp.


Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

[<u>IMAGES</u>]

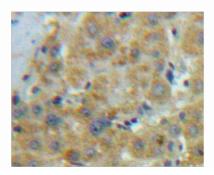

Used in DAB staining on fromalin fixed paraffin-embedded Liver tissue

[CONTENTS]

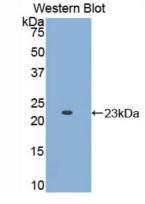

Used in Western Blot, Sample: Recombinant APOE, Human

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant APOE (Lys19~GIn181) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN $_3$ 0.02%.

[<u>STORAGE</u>]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.