

**PAA594Hu01****Polyclonal Antibody to Cytochrome C, Somatic (CYCS)****Organism Species: Homo sapiens (Human)*****Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

**[ PRODUCT INFORMATION ]****Immunogen:** CYCS, Human**Clonality:** Polyclonal**Host:** Rabbit**Immunoglobulin Type:** IgG**Purification:** Affinity Chromatography.**Applications:** WB, ICC, IHC-P, IHC-F, ELISA**Concentration:** 200µg/mL**UOM:** 100µg**[ IMMUNOGEN INFORMATION ]****Immunogen:** Recombinant CYCS (Glu5~Lys100) with N-terminal His-Tag expressed in *E.coli*.**Accession No.:** RPA594Hu01**[ ANTIBODY SPECIFICITY ]**

The antibody is a rabbit polyclonal antibody raised against CYCS. It has been selected for its ability to recognize CYCS in immunohistochemical staining and western blotting.

**[ APPLICATIONS ]**

Western blotting: 1:50-400

Immunocytochemistry in formalin fixed cells: 1:50-500

Immunohistochemistry in formalin fixed frozen section: 1:50-500

Immunohistochemistry in paraffin section: 1:10-100

Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

## [ CONTENTS ]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

## [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant CYCS (Glu5~Lys100) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.  
5uL per well when used in enhanced chemiluminescent (ECL).


**Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading Buffer:** 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN<sub>3</sub> 0.02%.


## [ STORAGE ]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.

## [ IMAGES ]

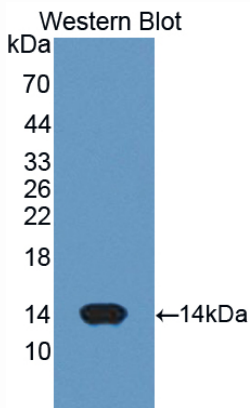

Used in Western Blot, Sample:  
Recombinant CYCS, Human

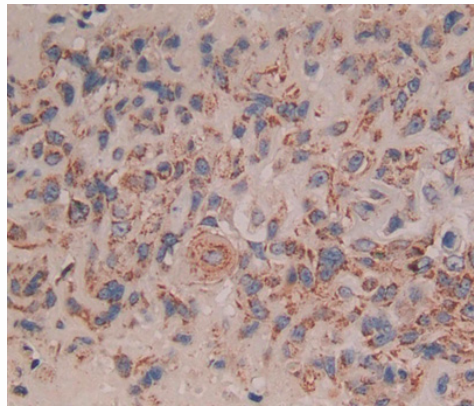

Used in DAB staining on fromalin fixed  
paraffin- embedded lung cancer tissue

**PAA594Hu01**

**Polyclonal Antibody to Cytochrome C, Somatic (CYCS)**

**Organism Species: Homo sapiens (Human)**


***Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY


NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---


---

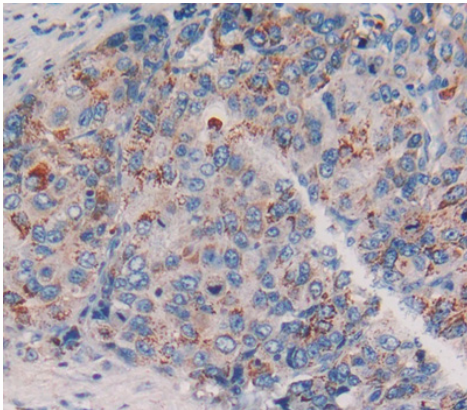


Used in DAB staining on formalin fixed  
paraffin- embedded liver cancer tissue

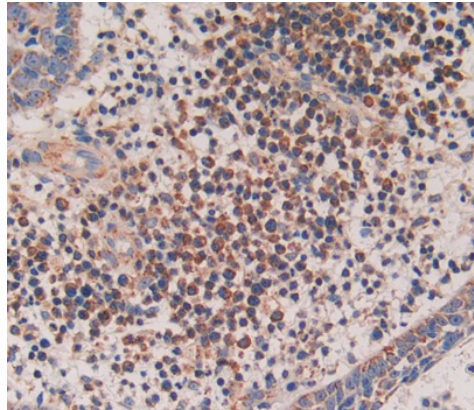

Used in DAB staining on formalin fixed  
paraffin- embedded skin cancer tissue

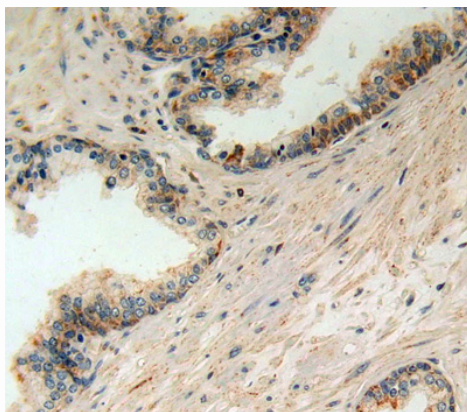

Used in DAB staining on formalin fixed  
paraffin- embedded prostate glands tissue

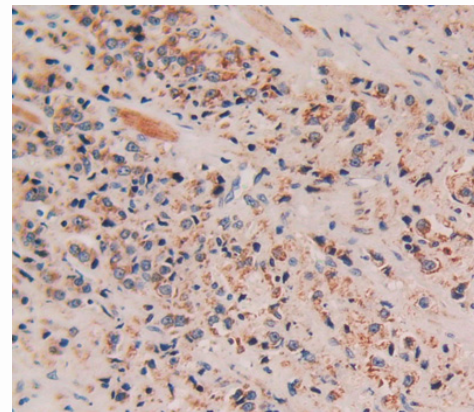

Used in DAB staining on formalin fixed paraffin-  
embedded prostate glands cancer tissue



**PAA594Hu01**

**Polyclonal Antibody to Cytochrome C, Somatic (CYCS)**

**Organism Species: Homo sapiens (Human)**


***Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY


NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---


---

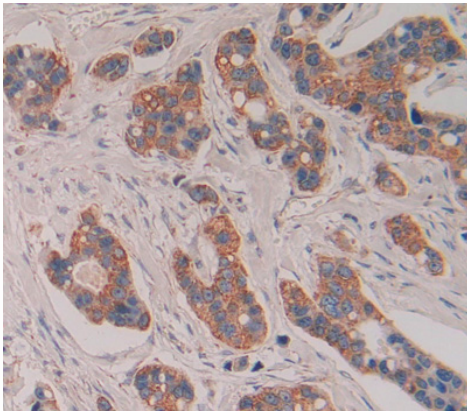


Used in DAB staining on formalin fixed  
paraffin- embedded breast cancer tissue

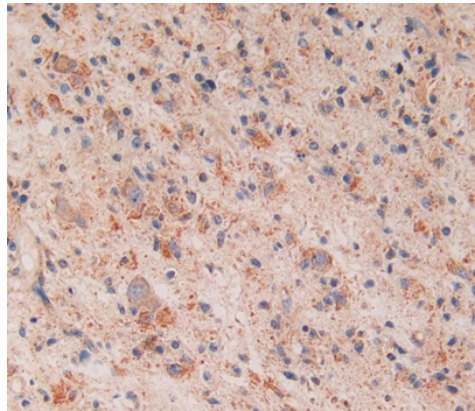

Used in DAB staining on formalin fixed  
paraffin- embedded glioma tissue

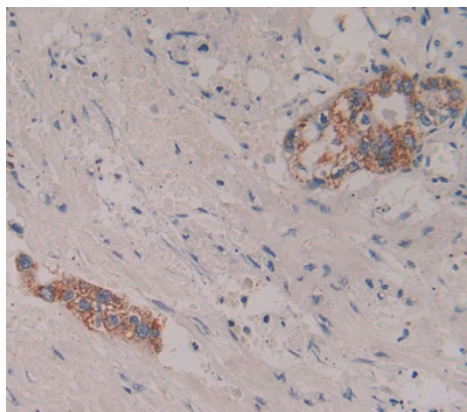

Used in DAB staining on formalin fixed  
paraffin- embedded pancreas cancer tissue

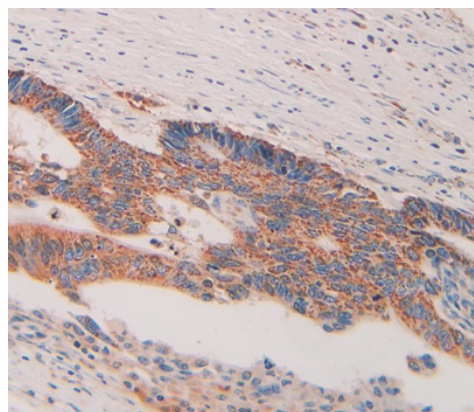

Used in DAB staining on formalin fixed  
paraffin- embedded rectum cancer tissue