

PAA624Hu01 Polyclonal Antibody to Hepatoma Derived Growth Factor (HDGF) Organism Species: Homo sapiens (Human) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

1304 Langham Creek Dr. Suite 226, Houston, TX 77084, USA | 001-888-960-7402 | www.cloud-cloue.us | mail@cloud-clone.as Export Processing Zone, Wuhan, Hubei 436056, PRC | 0086-000-880-0687 | www.cloud-clone.com | mail@cloud-clone.com

Cloud-Clone Corp.

[PROPERTIES]

Source: Polyclonal antibody preparation Host: Rabbit Purification: Antigen-specific Affinity Chromatography. Traits: Liquid Concentration: 200µg/mL UOM: 100µg Applications: WB; IHC; ICC; IP.

[IMMUNOGEN]

Immunogen: Recombinant HDGF (Tyr10~Leu240) expressed in *E.coli*. **Accession No.:** RPA624Hu01

[APPLICATIONS]

Western blotting: 0.5-5µg/mL Immunohistochemistry: 5-20µg/mL Immunocytochemistry: 5-20µg/mL Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant HDGF disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control.

Not used for other purposes.

Loading Buffer: 100mM Tris(pH6.8), 1% SDS, 150mM NaCl, 50% glycerol, 0.02% BPB, 50mM DTT and 0.02% NaN₃.

Coud-Clone Corp.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

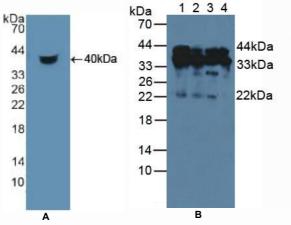


Figure 1. Western Blot
A. Sample: Recombinant HDGF, Human
B. Lane1: Human Hela Cells
Lane2: Human A549 Cells
Lane3: Human Jurkat Cells
Lane4: Human 293T Cells
Primary Ab: 3µg/mL Rabbit Anti-Human HDGF Ab
Second Ab: 1:2000 Dilution of HRP-Linked Guinea pig
Anti-Rabbit Ab (Catalog: SAA544Rb59)

Coud-Clone Corp.

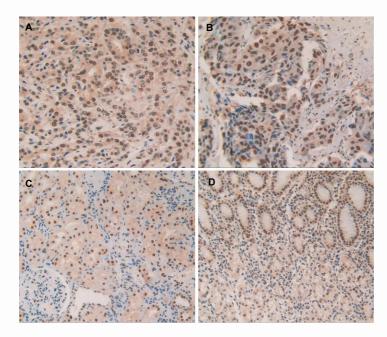


Figure 2. DAB staining on IHC-P

Samples:

- A. Human Prostate Gland Cancer Tissue
- B. Human Breast Cancer Tissue
- C. Human Kidney Tissue
- D. Human Stomach Tissue

1304 Langham Creek Dr. Soire 226, Houston, TX 77084, USA | 001-888-960-7402 | www.cloud-cloue.ox | mail@cloud-clone.ax Export Processing Zone, Wuhan, Hubei 430056, PRC | 0086-000-880-0687 | www.cloud-clone.com | mail@cloud-clone.com