PAB205Mu01 Polyclonal Antibody to Lymphocyte Cytosolic Protein 2 (LCP2) Organism Species: Mus musculus (Mouse) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

#### [PRODUCT INFORMATION]

Immunogen: LCP2, Mouse Clonality: Polyclonal Host: Rabbit Immunoglobulin Type: IgG Purification: Affinity Chromatography. Applications: WB, ICC, IHC-P, IHC-F, ELISA Concentration: 200µg/mL UOM: 100µg

#### [ IMMUNOGEN INFORMATION ]

Immunogen: Recombinant LCP2 (Thr202~Val463) expressed in E.coli.

Accession No.: RPB205Mu01

**Sequence:** The target protein is fused with two N-terminal Tags, His-tag and T7-tag and its sequence is listed below.

MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF- TGRNHSPLS PPHPNHEEPS RSGNNKTAKL PAPSIDRSTK PPLDRSLAPL DREPFILGKK PPFSDKPSAP LGREHLPKIQ KPPLPPAMDR HERNERLGPV TTRKPPVPRH GRGPDRREND EDDVHQRPLP QPSLPSMSSN TFPSRSVQPS SKNTFPLAHM PGAFSESNIG FQQSASLPPY FSQGPGNRPP LRSEGRNLPL PVPNRPQPPS PGEEETPLDE EWYVSYITRP EAEAALRKIN QDGTFLVRDS SKKTANNPYV LMV

#### [ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against LCP2. It has been selected for its ability to recognize LCP2 in immunohistochemical staining and western blotting.

#### [APPLICATIONS]


Western blotting: 1:100-400 Immunocytochemistry in formalin fixed cells: 1:100-500

Designed by Cloud-Clone Corp., Assembled by Uscn Life Science Inc. ISO9001:2008; ISO13485:2003 11271 Richmond Avenue Suite H104. Houston, TX 77082, USA | Toll free: 001-888-960-7402 | Fax: 001-832-538-0068 | Http://www.cloud-clone.us | E-mail: mail@cloud-clone.us Export Processing Zone Building F. Wuhan, Hubei 430056, PRC | Toll free: 0086-800-880-0667 | Fax: 0086-27-8425-9551 | Http://www.uscnk.com | E-mail: mail@cloud-clone.us

# Cloud-Clone Corp.

Immunohistochemistry in formalin fixed frozen section: 1:100-500 Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

#### [<u>IMAGES</u>]

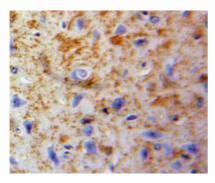

Used in DAB staining on fromalin fixed

paraffin-embedded Brain tissue

#### 

Recombinant LCP2, Mouse

### [ CONTENTS ]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

# [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant LCP2 (Thr202~Val463) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN $_3$  0.02%.

# [STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.

Designed by Cloud-Clone Corp., Assembled by Uscn Life Science Inc. ISO9001:2008; ISO13485:2003 11271 Richmond Avenue Suite H104, Houston, TX 77082, USA | Toll free: 001-888-960-7402 | Fax: 001-832-538-0088 | Http://www.cloud-clone.us | E-mail: mail@cloud-clone.us Export Processing Zone Building F. Wuhan, Hubel 430056, PRC | Toll free: 0086-800-880-0687 | Fax: 0086-27-8425-9551 | Http://www.uscnk.com | E-mail: mail@uscnk.com