

PAB198Ra01

Polyclonal Antibody to Nucleoporin 88kDa (NUP88)

Organism Species: Rattus norvegicus (Rat)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PRODUCT INFORMATION]

Immunogen: NUP88, Rat **Purification:** Affinity Chromatography.

Clonality: Polyclonal Applications: WB, ICC, IHC-P, IHC-F, ELISA

Host: Rabbit Concentration: 200µg/mL

Immunoglobulin Type: IgG **UOM**: 100μg

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant NUP88 (Thr159~Ile422) expressed in E.coli.

Accession No.: RPB198Ra01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and

T7-tag and its sequence is listed below.

MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF- TI PIAERFFTSS TSLTLKHAAW YPSEMLDPHI VLLTSDNVIR IYSLREPQTP TKVIVLSEAE EESLILNKGR AYTASLGETA VAFDFGPLVT VSKNMFEQKD REAVAYPLYI LYENGETFLT YVSLLHSPGN IGKLLGPLPM HPAAEDNYGY DACAILCLPC VPNILVIATE SGMLYHCVVL EGEEDDDQTL EKSWDPRADL IPSLYVFECV ELELALKLAS AEDDPFASDF SCPIKLHRDP KCPSRYHCSH EAGVHSVGLT WI

[ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against NUP88. It has been selected for its ability to recognize NUP88 in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500


Immunohistochemistry in formalin fixed frozen section: 1:100-500




Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

[IMAGES]

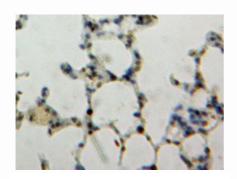

Used in DAB staining on fromalin fixed paraffin-embedded Lung tissue

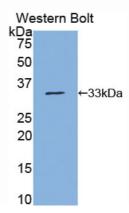

Used in Western Blot, Sample: Recombinant NUP88, Rat

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant NUP88 (Thr159~Ile422) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate. 5uL per well when used in enhanced chemilumescent (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN₃ 0.02%.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.