

PAB374Hu01

Polyclonal Antibody to Small Nuclear Ribonucleoprotein Polypeptide C (SNRPC)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

## [ PRODUCT INFORMATION ]

**Immunogen:** SNRPC, Human **Purification:** Affinity Chromatography.

Clonality: Polyclonal Applications: WB, ICC, IHC-P, IHC-F, ELISA

Host: Rabbit Concentration: 200µg/mL

**Immunoglobulin Type**: IgG **UOM**: 100μg

## [ IMMUNOGEN INFORMATION ]

Immunogen: Recombinant SNRPC (Pro2~Arg159) expressed in E.coli.

Accession No.: RPB374Hu01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and

T7-tag and its sequence is listed below.

MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF- PKFYCDYCD TYLTHDSPSV RKTHCSGRKH KENVKDYYQK WMEEQAQSLI DKTTAAFQQG KIPPTPFSAP PPAGAMIPPP PSLPGPPRPG MMPAPHMGGP PMMPMMGPPP PGMMPVGPAP GMRPPMGGHM

PMMPGPPMMR PPARPMMVPT RPGMTRPDR

## [ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against SNRPC. It has been selected for its ability to recognize SNRPC in immunohistochemical staining and western blotting.

## [APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500


Immunohistochemistry in paraffin section: 1:50-200




Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

#### [ IMAGES ]

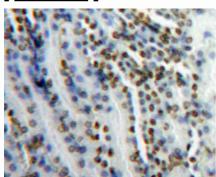

Used in DAB staining on fromalin fixed paraffin-embedded Kidney tissue

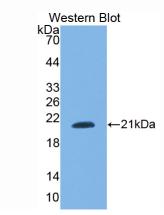

Used in Western Blot, Sample:
Recombinant SNRPC, Human

# [CONTENTS]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

# [QUALITY CONTROL]

**Content:** The quality control contains recombinant SNRPC (Pro2~Arg159) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL).

**Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading Buffer:** 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN<sub>3</sub> 0.02%.

## [STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.