PAB949Po01 Polyclonal Antibody to Transforming Growth Factor Beta 3 (TGFb3) Organism Species: Sus scrofa; Porcine (Pig) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PRODUCT INFORMATION]

C Cloud-Clone Corp.

Immunogen: TGFb3, Porcine Clonality: Polyclonal Host: Rabbit Immunoglobulin Type: IgG Purification: Affinity Chromatography. Applications: WB, ICC, IHC-P, IHC-F, ELISA Concentration: 200µg/mL UOM: 100µg

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant TGFb3 (Ala298~Ser409) with two N-terminal Tags, His-tag and T7-tag expressed in *E.coli*.

Accession No.: RPB949Po01

[ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against TGFb3. It has been selected for its ability to recognize TGFb3 in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:50-400 Immunocytochemistry in formalin fixed cells: 1:50-500 Immunohistochemistry in formalin fixed frozen section: 1:50-500 Immunohistochemistry in paraffin section: 1:10-100 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

Coud-Clone Corp.

[<u>CONTENTS</u>]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant TGFb3 (Ala298~Ser409) disposed in loading buffer.


Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN₃ 0.02%.

[<u>STORAGE</u>]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.

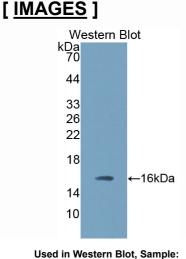

Recombinant TGFb3, Porcine

1304 Langham Creek Dr, Suite 226, Houston, TX 77084, USA | 001-888-960-7402 | www.cloud-clone.us | mail@cloud-clone.us Export Processing Zone, Wuhan, Hubei 430056, PRC | 0086-800-880-0687 | www.cloud-clone.com | mail@cloud-clone.com