

PAA213Bo01
Polyclonal Antibody
To Triggering Receptor Expressed On Myeloid Cells 1 (TREM1)
Organism Species: Bos taurus; Bovine (Cattle)
Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

[PRODUCT INFORMATION]

Immunogen: TREM1, Bovine

Clonality: Polyclonal

Host: Rabbit


Immunoglobulin Type: IgG

Purification: Affinity Chromatography.

Applications: WB, ICC, IHC-P, IHC-F, ELISA

Concentration: 200µg/mL

UOM: 100μg

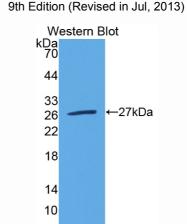

Sample: Recombinant TREM1, Bovine

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant TREM1 (Ala22~Ser232) expressed in E.coli.

Accession No.: RPA213Bo01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and

T7-tag and its sequence is listed below.

MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGS- AEVFEEKCT LAEGQTLKVS CPTNTNIYSN SQKAWQRLKD NGEVQTLAIT EGSSQVRVGK YFLEDIPSEG MLQIQMANLQ VEDSGLYRCV ILGPSDPIIL FHPVRLVVTK NSLGTPASDE YPCQVSVQNP TPLPVTTKLR PRPRPKPV TQPIPTSADR LSSPGFTVTP TNVTHVNRAP GISIIIPAAC GLLSKTLVFI GI FAVTHRSF AS



[ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against TREM1. It has been selected for its ability to recognize TREM1 in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant TREM1 (Ala22~Ser232) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL).

Note: The quality control is specifically manufactured as the positive control.

Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN₃ 0.02%.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.