

RPB874Ra02 10µg

Recombinant Coagulation Factor VII (F7)

Organism Species: Rattus norvegicus (Rat)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)



[PROPERTIES]

Source: Prokaryotic expression.

Host: E. coli

Residues: Ala315~Asp433
Tags: N-terminal His-Tag
Tissue Specificity: Plasma.

Subcellular Location: Secreted.

Purity: >98%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA,

1mM DTT, 0.01% sarcosyl, 5%Trehalose and Proclin300.

Original Concentration: 200ug/mL

Applications: SDS-PAGE; WB; ELISA; IP; CoIP; Purification; Amine Reactive

Labeling.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.4

Predicted Molecular Mass: 14.4kDa

Accurate Molecular Mass: 14kDa as determined by SDS-PAGE reducing conditions.

[<u>USAGE</u>]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

ASIRFS RVSGWGQLLD RGATALELMV IEVPRLMTQD CLEHAKHSAN TPRITENMFC AGYMDGTKDA CKGDSGGPHA THYHGTWYLT GVVSWGEGCA AIGHIGVYTR VSQYIDWLVK YMD

[IDENTIFICATION]

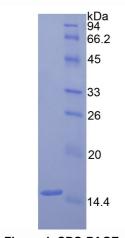


Figure 1. SDS-PAGE