

RPD696Hu01 500µg
Recombinant Diazepam Binding Inhibitor (DBI)
Organism Species: *Homo sapiens (Human)*
Instruction manual

FOR RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Met1~Ile104

Tags: N-terminal His Tag

Subcellular Location: Golgi apparatus

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% Sarcosyl, 5% Trehalose.

Original Concentration: 300µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.4

Predicted Molecular Mass: 15.5kDa

Accurate Molecular Mass: 16kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

MWGDLWLLPP ASANPGTGTE AEFKAAEEV RHLKTKPSDE EMLFIYGHYK
 QATVGDINTE RPGMLDFTGK AKWDAWNEK GTSKEDAMKA YINKVEELKK
 KYGI

[IDENTIFICATION]

ATGTGGGGGACCTCTGGCTCTCTGGCTCTCTGCAATCCGAGCCTGGGACAGAGGCTGAGTTTGGGAACTGAGAGGAGSTTAGCCOCTTAAGACCAACCACTGATGAGGAGTGTGTTTCACTATGACCACTTCAAAACAAGCAACTGTGGGACAAATAAGAGAACGGCCCGCATGTTGGACTCAAGGCT
 H V G D I V L L P P A S A N P G T G T E A E F E K A A E E V R H L K T K P S D E E M L F I Y G H Y K Q A T V G D I N T E R P G H L D F T G

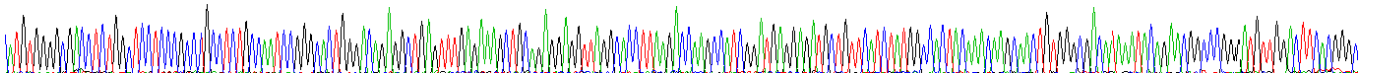


Figure. Gene Sequencing (Extract)

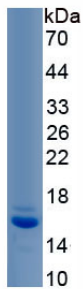


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.