

RPD017Ra01 10 μ g
Recombinant Follicle Stimulating Hormone Beta (FSH β)
Organism Species: *Rattus norvegicus* (Rat)
Instruction manual

FOR RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[**PROPERTIES**]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Ser21~Glu130

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.3

Predicted Molecular Mass: 16.2kDa

Accurate Molecular Mass: 16kDa as determined by SDS-PAGE reducing conditions.

[**USAGE**]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[**STORAGE AND STABILITY**]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[**SEQUENCE**]

SCELTNITIS VEKEECRFCI SINTTWCEGY
CYTRDLVYKD PARPNTQKVC TFKELVYETI RLPGCARHSD SLYTYPVATE
CHCGKCDSDS TDCTVRGLGP SYCSFGEMKE

[IDENTIFICATION]

TCGAGCGATGATGACGCGAGCGCTTCAGACGATTCGGGGATCATTGGTCCGAGCGCTGCAAGATCGAGGCTTACCTGGACCTTACCTGGCGGATCGCTTCACGATGCGGATGAGCGGACGATCGATGAGCGCTGCGAT
SCELTNITISVEKEECRFCISINTTWCEGYCYTRDLVYKDPARPNTQKVCTFKELVYETIRLPGCARHSDSLYTYPVATECHCGKCDSDSDCTVRGLGPSY


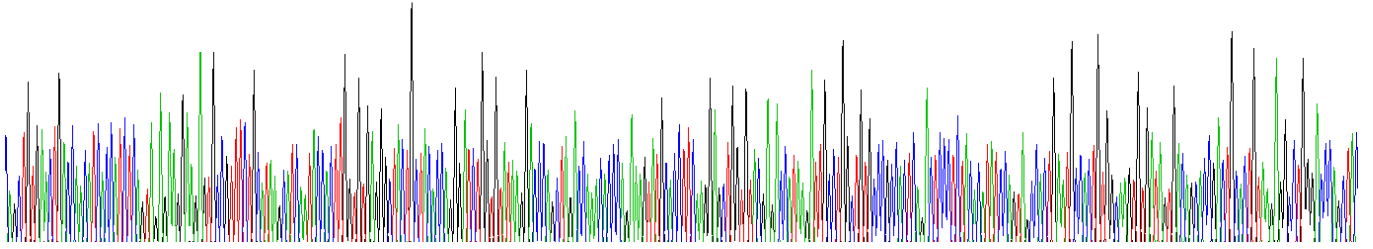


Figure. Gene Sequencing (Extract)

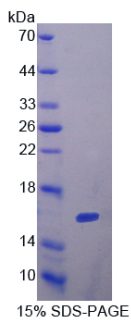

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.