

RPE531Hu01 50µg Recombinant Gasdermin A (GSDMA) Organism Species: *Homo sapiens (Human) Instruction manual*

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

Cloud-Clone Corp.

[PROPERTIES]

Source: Prokaryotic expression Host: *E.coli* Residues: Met1~Val336 Tags: N-terminal His Tag Subcellular Location: Cytoplasm Purity: > 90% Traits: Freeze-dried powder Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose . Original Concentration: 200µg/mL Applications: Positive Control; Immunogen; SDS-PAGE; WB. (May be suitable for use in other assays to be determined by the end user.) Predicted isoelectric point: 5.6 Predicted Molecular Mass: 41.1kDa

Accurate Molecular Mass: 43kDa as determined by SDS-PAGE reducing conditions.

[<u>USAGE</u>]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

Cond-Clone Corp.

MTMFENVTRA LARQLNPRGD LTPLDSLIDF KRFHPFCLVL RKRKSTLFWG ARYVRTDYT LLDVLEPGSS PSDPTDTGNF GFKNMLDTRV EGDVDVPKTV KVKGTAGLS QNSTLEVQTL SVAPKALETV QERKLAADHP FLKEMQDQGE NLYVVMEVV ETVQEVTLER AGKAEACFSL PFFAPLGLQG SINHKEAVTI PKGCVLAFR VRQLMVKGKD EWDIPHICND NMQTFPPGEK SGEEKVILIQ ASDVGDVHE GFRTLKEEVQ RETQQVEKLS RVGQSSLLSS LSKLLGKKKE LQDLELALE GALDKGHEVT LEALPKDVLL SKEAVGAILY FV

[IDENTIFICATION]

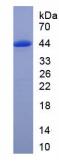


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.