

RPJ822Hu01 1

Recombinant Inhibitor Of Kappa-Light Polypeptide Gene Enhancer In B-Cells Kinase

Beta (IkBKb)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

# Cond-Clone Corp.

## [PROPERTIES]

Source: Prokaryotic expression

Host: E.coli

Residues: Asp246~Val528

Tags: N-terminal His and GST Tag

Subcellular Location: Nucleus, Cytoplasm

**Purity:** > 90%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA, 1mM DTT,

0.01% SKL, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.4

Predicted Molecular Mass: 62.3kDa

Accurate Molecular Mass: 62kDa as determined by SDS-PAGE reducing conditions.

#### [<u>USAGE</u>]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

#### [ STORAGE AND STABILITY ]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.


**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.



### [<u>SEQUENCE</u>]

DLNGT VKFSSSLPYP NNLNSVLAER LEKWLQLMLM WHPRQRGTDP TYGPNGCFKA LDDILNLKLV HILNMVTGTI HTYPVTEDES LQSLKARIQQ DTGIPEEDQE LLQEAGLALI PDKPATQCIS DGKLNEGHTL DMDLVFLFDN SKITYETQIS PRPQPESVSC ILQEPKRNLA FFQLRKVWGQ VWHSIQTLKE DCNRLQQGQR AAMMNLLRNN SCLSKMKNSM ASMSQQLKAK LDFFKTSIQI DLEKYSEQTE FGITSDKLLL AWREMEQAVE LCGRENEV

## [IDENTIFICATION]

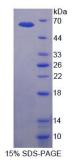

#### [IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.