

RPC180Mu01 1

Recombinant Interferon Regulatory Factor 2 (IRF2)

Organism Species: *Mus musculus* (Mouse)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[**PROPERTIES**]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Met1~Cys349

Tags: N-terminal His Tag

Subcellular Location: Nucleus

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 100µg

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.2

Predicted Molecular Mass: 43.1kDa

Accurate Molecular Mass: 55kDa as determined by SDS-PAGE reducing conditions.

Phenomenon explanation:

The possible reasons that the actual band size differs from the predicted are as follows:

1. Splice variants: Alternative splicing may create different sized proteins from the same gene.
2. Relative charge: The composition of amino acids may affect the charge of the protein.
3. Post-translational modification: Phosphorylation, glycosylation, methylation etc.
4. Post-translation cleavage: Many proteins are synthesized as pro-proteins, and then cleaved to give the active form.
5. Polymerization of the target protein: Dimerization, multimerization etc.

[**USAGE**]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[**STORAGE AND STABILITY**]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

```

MGSSHHHHH  SSSLVPRGSH  MASMTGGQQM  GRGS-MPVERMRMP  WLEEQINSNT
IPGLKWLNKE  KKIFQIPWMH  AARHGWDVEK  DAPLFRNWAI
HTGKHQPGID  KPDPKTWKAN  FRCAMNSLPD  IEEVKDRSIK  KGNNAFRVYR  MLPLSERPSK
KGKKPKTEKE  ERVKHIKQEP  VESSLGLSNG  VSGFSPEYAV  L TSAIKNEVD  STVNIIVVGQ
SHLDSNIEDQ  EIVTNPPDIC  QVVEVTTESD  DQPVSMSELY  PLQISPVSSY  AESETTDSVA
SDEENAEGRP  HWRKRSIEGK  QYLSNMGTRN  TYLLPSMATF  VTSNKPDQV  TIKEDSCPMP
YNSSWPPFTD  LPLPAPVTPT  PSSSRPDRET  RASVIKKTSD  ITQARVKSC
    
```

[IDENTIFICATION]


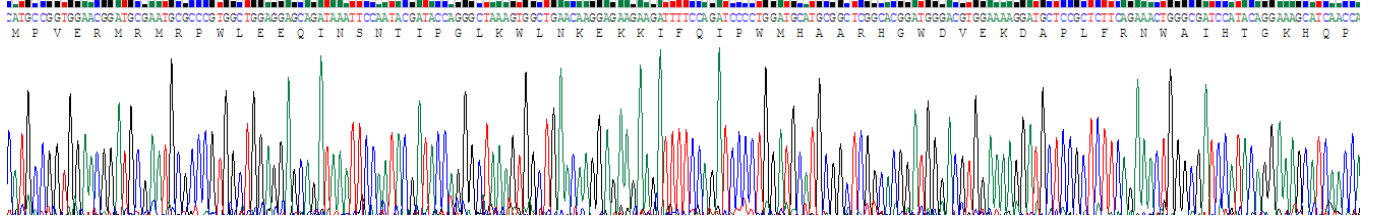


Figure. Gene Sequencing (Extract)

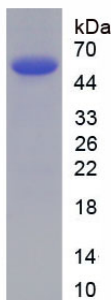

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.