

RPB630Hu01 50µg

Recombinant Latent Transforming Growth Factor Beta Binding Protein 2 (LTBP2)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Ser522~Thr814

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 95%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA, 1mM DTT, 0.01% SKL, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.0

Predicted Molecular Mass: 35.0kDa

Accurate Molecular Mass: 38kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[**SEQUENCE**]

		SLWDSNNIP	ARSGEPPRPL	PPAAPRRPRL
LGRCYLNTVN	GQCANPLLEL	TTQEDCCGSV	GAFWGVTLCA	PCPPRPASPV
IENGQLECPQ	GYKRLNLTHC	QDINECLTLG	LCKDAECVNT	RGSYLCTCRP
GLMLDPSRSR	CVSDKAISML	QGLCYRSLGP	GTCTLPLAQR	ITKQICCCSR
VGKAWGSECE	KCPLPGTEAF	REICPAGHGY	TYASSDIRLS	MRKAEEEELA
RPPREQGQRS	SGALPGPAER	QPLRVVTDTW	LEAGTIPDKG	DSQAGQVTTS
VTHAPAWVTG	NATT			

[**IDENTIFICATION**]

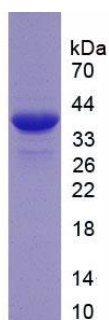


Figure. SDS-PAGE

[**IMPORTANT NOTE**]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.