

RPC922Hu01 100µg

Recombinant Platelet Derived Growth Factor Receptor Like Protein (PDGFRL)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Gln22~Ser375

Tags: Two N-terminal Tags, His-tag and GST-tag

Accession: Q15198

Host: E. coli

Subcellular Location: Secreted.

Purity: >95%

Endotoxin Level: <1.0EU per 1μg (determined by the LAL method).

Formulation: Supplied as lyophilized form in PBS,


pH7.4, containing 5% trehalose, 0.01% sarcosyl.

Predicted isoelectric point: 9.0

Predicted Molecular Mass: 69.5kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

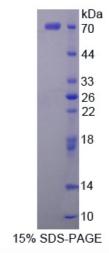

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.



[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

QHLPKNKRP KEPGENRIKP TNKKVKPKIP KMKDRDSANS APKTQSIMMQ VLDKGRFQKP AATLSLLAGQ TVELRCKGSR IGWSYPAYLD TFKDSRLSVK QNERYGQLTL VNSTSADTGE FSCWVQLCSG YICRKDEAKT GSTYIFFTEK GELFVPSPSY FDVVYLNPDR QAVVPCRVTV LSAKVTLHRE FPAKEIPANG TDIVYDMKRG FVYLQPHSEH QGVVYCRAEA GGRSQISVKY QLLYVAVPSG PPSTTILASS NKVKSGDDIS VLCTVLGEPD VEVEFTWIFP GQKDERPVTI QDTWRLIHRG LGHTTRISQS VITVEDFETI DAGYYICTAQ NLQGQTTVAT TVEFS