RPH871Mu01 100µg **Recombinant Pleckstrin Homology Domain** Containing Family A, Member 1 (PLEKHA1) **Organism Species: Mus musculus (Mouse)** Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Thr105~Leu351 Tags: Two N-terminal Tags, His-tag and GST-tag Accession: Q8BUL6 Host: E. coli **Subcellular Location:** Cytoplasm. Cell membrane; Peripheral membrane protein. Nucleus. **Purity:** >90% Endotoxin Level: <1.0EU per 1µg (determined by the LAL method). Formulation: Supplied as lyophilized form in PBS, pH7.4, containing 5% trehalose, 0.01% sarcosyl. Predicted isoelectric point: 8.8


Predicted Molecular Mass: 57.3kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.

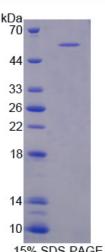

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[<u>SEQUENCES</u>]

The sequence of the target protein is listed below.

TPTQKE EVNECGESLD RNNLKRSQSH LPYFAPKPPS DSAVIKAGYC VKQGAVMKNW KRRYFQLDEN TIGYFKSELE KEPLRVIPLK EVHKVQECKQ SDIMMRDNLF EIVTTSRTFY VQADSPEEMH SWIKAVSGAI VAQRGPGRSS SSEHSGWPSE PKFTVRPAGA AATASHSPAS RSHSLVSGCA MEKRGVYQSL AKVKPGNFKV QTVSPREPAS KVTEQALLRP QSKNGPQGED GEPVDLDDAS L