

RPG078Mu01 100ug

Recombinant Purinergic Receptor P2X, Ligand Gated Ion Channel 1 (P2RX1)

Organism Species: *Mus musculus* (Mouse)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Thr57~Gly338

Tags: N-terminal His Tag

Subcellular Location: Membrane

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: pH7.4, containing 0.01% SKL, 5% Trehalose

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.0

Predicted Molecular Mass: 35kDa

Accurate Molecular Mass: 36kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

TSSG LISSVSVKLLK GLAVTQLQGL GPQVWDVADY VFPAHGDSSF VVMTNFIMTP
QQAQGHCAEN PEGGICQDDS GCTPGKAERK AQGIRTGNCV PFNGTVKTCE
IFGWCPVEVD DKIPSPALLH EAENFTLFIK NSISFPRFKV NRRNLVEEVN
GTYMKKCLYH KILHPLCPVF SLGYVVRESG QDFRSLAEKG GVVGITIDWE
CDLDWHVRHC KPIYQFHGLY GEKNLSPGFN FRFARHFVQN GTNRRHLFKV
FGIRFDILVD GKAGKFDIIP TMTTIGSG

[IDENTIFICATION]

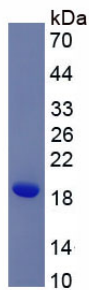

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.