

**RPE971Hu01 200µg**

**Recombinant Ribosomal Protein, Large, P0 (RPLP0)**

**Organism Species: *Homo sapiens* (Human)**

***Instruction manual***

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

13th Edition (Revised in Aug, 2023)

**[ PROPERTIES ]**

**Source:** Prokaryotic expression

**Host:** *E.coli*

**Residues:** Met1~Asp317

**Tags:** N-terminal His Tag

**Subcellular Location:** Nucleus, Cytoplasm

**Purity:** > 97%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 5% Trehalose .

**Original Concentration:** 500µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.7

**Predicted Molecular Mass:** 38.0kDa

**Accurate Molecular Mass:** 38kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]**

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

**[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ SEQUENCE ]**

MPREDRATWK SNYFLKIIQL LDDYPKCFIV GADNVGSKQM QQIRMSLRGK  
AVVLMGKNTM MRKAIRGHLE NNPALKLLP HIRGNVGFVF TKEDLTEIRD  
MLLANKVPAA ARAGAIAPCE VTVPAQNTGL GPEKTSFFQA LGITTKISRG  
TIEILSDVQL IKTGDKVGAS EATLLNMLNI SPFSFGLVIQ QVFDNGSIYN  
PEVLDITEET LHSRFLEGVR NVASVCLQIG YPTVASVPHS IINGYKRVLA  
LSVETDYTFP LAEKVKAFLA DPSAFVAAAP VAAATTAAPA AAAAPAKVEA  
KEESESEDED MGFGLFD

**[ IDENTIFICATION ]**

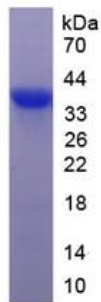


Figure. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.