

#### RPP155Bo01 100µg

Recombinant Wingless Type MMTV Integration Site Family, Member 3A (WNT3A)

Organism Species: Bos taurus; Bovine (Cattle)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

# [PROPERTIES]

Residues: Val19~Lys355

Tags: Two N-terminal Tags, His-tag and T7-tag

Accession: NM\_001206024.1

Host: *E. coli*Purity: >95%

**Endotoxin Level:** <1.0EU per 1µg (determined by the LAL method).


**Formulation:** Supplied as lyophilized form in PBS, pH7.4, containing 5% trehalose, 0.01% sarcosyl.

Predicted isoelectric point: 7.5

Predicted Molecular Mass: 41.4kDa

**Applications:** SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

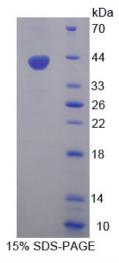

### [USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.



## [STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

## [SEQUENCES]

The sequence of the target protein is listed below.

VL AGYPIWWSLA LGQQYTSLGS QPLLCGSIPG LVPKQLRFCR NYIEIMPSVA EGVKLGIQEC QHQFRGRRWN CTTIDDSLAI FGPVLDKATR ESAFVHAIAS AGVAFAVTRS CAEGTSTICG CDSHHKGPPG EGWKWGGCSE DADFGVLVSR EFADARENRPDARSAMNKHN NEAGRTTILD HMHLKCKCHG LSGSCEVKTC WWAQPDFRAI GDFLKDKYDS ASEMVVEKHR ESRGWVETLR AKYALFKPPT ERDLVYYENS PNFCEPNPET GSFGTRDRTC NVTSHGIDGC DLLCCGRGHN TRTEKRKEKC HCIFHWCCYV SCQECIRIYD VHTCK