### RPH468Hu01 100µg Recombinant p21 Protein Activated Kinase 2 (PAK2) Organism Species: Homo sapiens (Human) *Instruction manual*

#### FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

# [PROPERTIES]

Residues: Glu177~Ala419 (Accession # Q13177),with two N-terminal Tags, His-tag and T7-tag.Host: *E. coli*Subcellular Location: Cytoplasm. Nucleus.Cytoplasm, perinuclear region. Membrane;Lipid-anchor.Purity: >95%Endotoxin Level: <1.0EU per 1µg</td>(determined by the LAL method).Formulation: Supplied as lyophilized form in PBS,pH7.4, containing 5% trehalose, 0.01% sarcosyl.Predicted isoelectric point: 6.6

Predicted Molecular Mass: 30.7kDa


Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

### [ <u>USAGE</u> ]

Reconstitute in sterile PBS, pH7.2-pH7.4.

Designed by Cloud-Clone Corp., Assembled by Uscn Life Science Inc. ISO9001:2008; ISO13485:2003; C € 11271 Richmond Avenue Suite H104, Houston, TX 77082, USA | Toll free: 001-888-960-7402 | Fax: 001-832-538-0088 | Http://www.cloud-clone.us | E-mail: mail@cloud-clone.us Expant Processing Zane Building F. Wuhan, Hubei 430056, PRC | Toll free: 0086-800-880-0687 | Fax: 0086-27-8425-9551 | Http://www.uscnk.com | E-mail: mail@uscnk.com

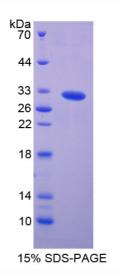

🗂 <u>Cloud-Clone</u> Corp.

# [ STORAGE AND STABILITY ]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

# [<u>SEQUENCES</u>]

The target protein is fused with two N-terminal Tags, His-tag and T7-tag, its sequence is listed below.

MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF- ETAP PVIAPRPDHT KSIYTRSVID PVPAPVGDSH VDGAAKSLDK QKKKTKMTDE EIMEKLRTIV SIGDPKKKYT RYEKIGQGAS GTVFTATDVA LGQEVAIKQI NLQKQPKKEL IINEILVMKE LKNPNIVNFL DSYLVGDELF VVMEYLAGGS LTDVVTETCM DEAQIAAVCR ECLQALEFLH ANQVIHRDIK SDNVLLGMEG SVKLTDFGFC AQITPEQSKR STMVGTPYWM APEVVTRKA