

MAB703Hu26

Monoclonal Antibody to TATA Binding Protein (TBP)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Monoclonal antibody preparation

Host: Mouse

Antibody isotype: IgG2a Kappa

Purification: Protein A + Protein G affinity chromatography

Clone number: C10

Traits: Liquid

Concentration: 1mg/mL

UOM: 100μL

Cross Reactivity: Mouse;Rat;Cavia;Rabbit;Porcine;Ovine;Gallus

Applications: WB

[IMMUNOGEN]

Immunogen: Recombinant TBP (Ile25~Thr339) expressed in *E.coli*

Accession No.: RPB703Hu01

[APPLICATIONS]

Western blotting: 0.01-2μg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.




Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the

expiration date under appropriate storage condition.

[IDENTIFICATION]

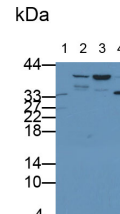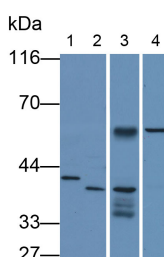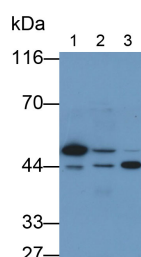

Western Blot; Samples: Lane1: Porcine Liver lysate; Lane2: HeLa cell lysate; Lane3: 293T cell lysate; Primary Ab: 1µg/ml Mouse Anti-Human TBP Antibody Second Ab: 0.2?g/ml HRP-Linked Caprine Anti-Mouse IgG Polyclonal Antibody (Catalog: SAA544Mu19)	Western Blot; Samples: Lane1: 293T cell lysate; Lane2: Cavia Testis lysate; Lane3: Mouse Testis lysate; Lane4: Gallus Kidney lysate; Primary Ab: 1µg/ml Mouse Anti-Human TBP Antibody Second Ab: 0.2?g/ml HRP-Linked Caprine Anti-Mouse IgG Polyclonal Antibody (Catalog: SAA544Mu19)	Western Blot; Samples: Lane1: Rabbit Kidney lysate; Lane2: Ovine Testis lysate; Lane3: Porcine Kidney lysate; Lane4: Rat Testis lysate; Primary Ab: 1µg/ml Mouse Anti-Human TBP Antibody Second Ab: 0.2?g/ml HRP-Linked Caprine Anti-Mouse IgG Polyclonal Antibody (Catalog: SAA544Mu19)
--	---	--

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.