

**RPB406Mu02 10 $\mu$ g**  
**Recombinant Lipopolysaccharide Binding Protein (LBP)**  
**Organism Species: *Mus musculus* (Mouse)**  
***Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

10th Edition (Revised in Jan, 2014)

## **[ PROPERTIES ]**

**Residues:** Gly25~Val481

**Tags:** Two N-terminal Tags, His-tag and T7-tag

**Accession:** Q61805

**Host:** *E. coli*

**Subcellular Location:** Secreted.

**Purity:** >90%

**Endotoxin Level:** <1.0EU per 1 $\mu$ g

(determined by the LAL method).


**Formulation:** Supplied as lyophilized form in 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA, 1mM DTT, 0.01% sarcosyl, 5% trehalose, and preservative.

**Predicted isoelectric point:** 8.7

**Predicted Molecular Mass:** 54.5kDa

**Applications:** SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

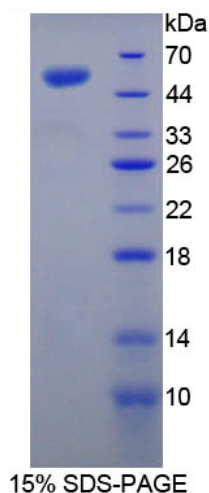

## **[ USAGE ]**

Reconstitute in sterile ddH<sub>2</sub>O.

## **[ STORAGE AND STABILITY ]**

**Storage: Avoid repeated freeze/thaw cycles.**

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

## **[ SEQUENCES ]**

The sequence of the target protein is listed below.

GVNPGV VARITDKGLA YAAKEGLVAL QRELYKITLP DFSGDFKIKV VGRGQYEFHS  
LEIQNCELRG SSLKLLPGQG LSLAISDSSI GVRGKWKVRK SFLKLHGSFD LDVKGVTISV  
DLLLGMDPSG RPTVSASGCS SRICDLVDHI SGNVGWLLNL FHNQIESKLQ KVLENKVCES  
IQKSVTSDLQ PYLQTLPVTA EIDNVLGIDY SLVAAPQAKA QVLDVMFKGE IFNRNHRSPV  
ATPTPTMSLP EDSKQMVYFA ISDYAFNIAS RVYHQAGYLN FSITDDMLPH DSGIRLNTKA  
FRPFTPQIYK KYPDMKLELL GTVVSAPILN VSPGNLSLAP QMEIEGFVIL PTSAREPVFR  
LGVVTNVFAS LTFNNSKVTG MLHPDKAQVR LIESKVGMFN VNLQFAFLNY YLLNSLYPDV  
NAELAQQGFPL PLPRHIQLHD LDFQIRKDFL YLGANVQYMR V

## **[ REFERENCES ]**

1. Lengacher S., *et al.* (1995) J. Inflamm. 47:165-172.
2. Eckert J.K., *et al.* (2013) Immunity 39:647-660.
3. Heinrich J.M., *et al.* (2001) J. Immunol. 167:1624-1628.
4. Uesugi T., *et al.* (2002) J. Immunol. 168:2963-2969.