

**RPJ828Hu01 20µg**

**Recombinant Phosphoinositide-3-Kinase Class 3 (PIK3C3)**

**Organism Species: *Homo sapiens (Human)***

***Instruction manual***

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

13th Edition (Revised in Aug, 2023)

**[ PROPERTIES ]**

**Source:** Prokaryotic expression

**Host:** *E.coli*

**Residues:** Tyr631~Trp885

**Tags:** N-terminal His Tag

**Subcellular Location:** Plasma, Lysosome

**Purity:** > 90%

**Traits:** Freeze-dried powder

**Buffer formulation:** 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA, 1mM DTT, 0.01% Sarcosyl, 5% Trehalose .

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.6

**Predicted Molecular Mass:** 32.9kDa

**Accurate Molecular Mass:** 33kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]**

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

**[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ SEQUENCE ]**

```
YPVIFKHGDD LRQDQLILQI
ISLMDKLLRK ENLDLKLTPY KVLATSTKHG FMQFIQSVPV AEVLDTEGSI
QNFFRKYAPS ENGPNGISAE VMDTYVKSCA GYCVITYILG VGDRHLDNLL
LTKTGKLFHI DFGYILGRDP KPLPPPMKLN KEMVEGMGGT QSEQYQEFRK
QCYTAFHLR RYSNLILNLF SLMVDANIPD IALEPDKTVK KVQDKFRLDL
SDEEAVHYMQ SLIDESVHAL FAAVVEQIHK FAQYW
```

**[ IDENTIFICATION ]**

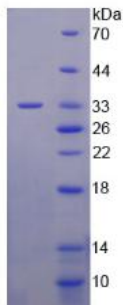


Figure. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.